


Please note:


- ❗ For this test, you'll need a standard adhesive strip (e.g. scotch tape).
- ❗ The adhesive strip must be transparent.
- ❗ Avoid contamination with fecal matter. This will affect the assessment.



- ➔ Remove the specimen slide from the protective tube.

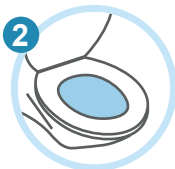

- ➔ First thing in the morning before going to the toilet, apply a transparent sticky strip with the sticky side up and press against your anus.

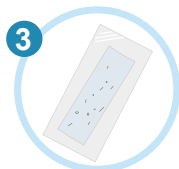


- ➔ Then, stick the adhesive tape onto the slide and put this back into the protective tube, closing it tightly



- ➔ Place the tube into the clear protective bag and seal it.
- ➔ Complete the accompanying sample submission form and sign it.

Your kit contains:

- ✓ These instructions
- 1 specimen slide in a protective tube
- 1 sample submission form
- 1 transparent bag
- 1 return envelop

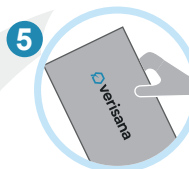

- ➔ Place the bag and the sample form in the grey return envelope and post it back to our laboratory.